



MELISSA BASIC - *Mattress overlay*

An affordable mattress overlay with pressure-relieving properties

The Melissa Basic Mattress overlay is a affordable mattress overlay designed to reduce the risk of pressure ulcers up to and including category 1. A 6 cm overlay mattress with pressure relieving properties, this mattress is a good complement to your existing homecare mattress. In addition, with the ability to flip and rotate the core, you get long-term use while maintaining quality.

- **HIGH RESILIENT COLD FOAM.** The Melissa Basic mattress overlay is made from high resilient cold foam which provides a durable and stable mattress.
- **REVERSIBLE FOAM CORE.** The foam core of this mattress overlay can be both flipped and rotated without altering its properties. By doing so, the lifespan of your Melissa 12 mattress overlay will be extended.
- **WATERPROOF COVER WITH ELASTIC.** The mattress overlay comes with a noiseless RX cover. The cover is supplied with elastic bands at the corners to keep the mattress overlay in place.




MELISSA BASIC - *Mattress overlay*

A GREAT COMPLEMENT TO YOUR HEALTHCARE MATTRESS.

Melissa Basic mattress overlay is designed to provide pressure relief and comfort and can also be used as a base and spacer mattress. The two-stretch hygienic cover is also designed to reduce the risk of shearing.

Waterproof cover with good properties

Melissa Basic mattress overlay comes with a cover from the RX series. A series of covers made of 2-stretch polyurethane fabric that is both breathable and noiseless. Depending on the choice of cover, the cover has zippers on three or four sides. All covers have a drip shield covering the zipper to facilitate hygiene.

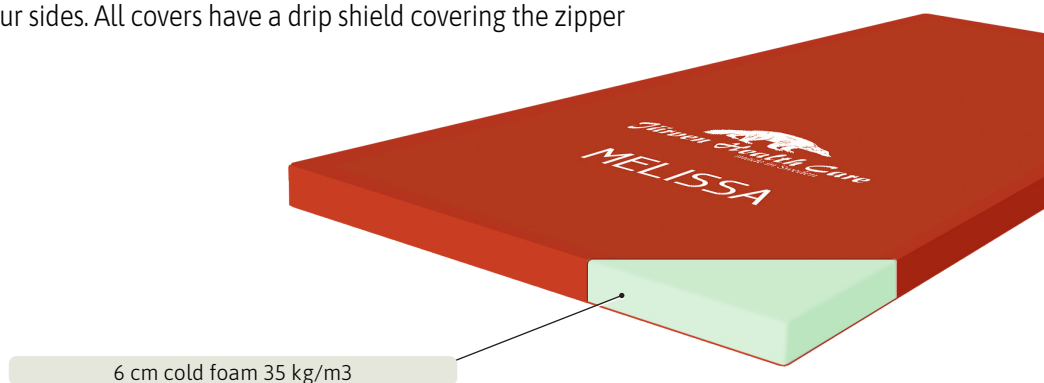

STANDARD FOR THE OBJECTIVE EVALUATION OF THE PRESSURE-RELIEVING PROPERTIES OF MATTRESSES.

In order to obtain an objective assessment of the pressure relieving properties of mattresses, RISE has developed a standard where the actual counterpressure of a mattress is measured by means of an indentational body. The measurement method is described in the Swedish standard SS 876 00 13.

Melissa Basic mattress overlay has been tested by RISE* according to method SS 876 00 13 and is pressure relieving up to wound category 2.

The test results are presented below.

| Load capacity | Extra low | Low | Medium | High |
|----------------------|-----------|-------------|-------------|-------------|
| Large indenter 300 N | < 30,0 | 30,0 - 34,9 | 35,0 - 39,9 | 40,0 - 44,9 |
| Large indenter 600 N | < 45,0 | 45,0 - 52,9 | 53,0 - 60,9 | 61,0 - 68,9 |
| Small indenter 40 N | < 45,0 | 45,0 - 52,9 | 53,0 - 60,9 | 61,0 - 68,9 |

| Mattress | Large indenter 300 N | Large indenter 600 N | Small indenter 40 N |
|---------------------------------|----------------------|----------------------|---------------------|
| Melissa Basic, mattress overlay | 34 | 52 | 60 |

*RISE (Research Institutes of Sweden), a Swedish independent research institute - www.ri.se

PRODUCT SPECIFICATIONS



| | | | |
|------------------------|--|--------------------|---|
| Materials | 6 cm cold foam 50 kg/m ³ | User weight | Up to and including 180 kg. |
| Sizes | 200/205/210 x 80/85/90/105/120 x 6 cm (For other sizes contact Järven Health Care.) | Fire safety | SS 87600 01, normal risk (test method EN 597-1 & 2) |
| Wound category | Up to and including category 1. | Recycling | Mattress with cover is fully recyclable. Sort as combustible waste. Disassembly is not required. |
| Mattress weight | approx. 4 kg depending on length and width. | Marking | Year of manufacture, month, day, EAN, care instructions, contact information, standards met and bed location. |



Bromsvägen 3 | SE-891 60 Örnsköldsvik | SWEDEN
Tel: +46 (0)660 - 22 28 50 | Fax: +46 (0)660 - 824 10
www.jarven.se | healthcare@jarven.se